

Capitainer® B10 Vanadate

50 pcs box.

Intended purpose:

Capitainer®B10 Vanadate sampling cards are used to collect and save a specific volume of capillary blood as a dried blood stain.

Material:

Sample disc:

Ahlstrom grade 222 Including
50 µg NaVO₃

Sampling card:

PET, acrylic adhesives, paper,
cardboard

Dimension:

- Sampling card: 43 x 84 mm
- Primary package: 80 x 120 mm
- Box 50 bundle: 85 x 125 x 230 mm

Product performance:

The product generates 2 x 10µl of dried whole blood sample from two drops of blood with a precision of < 2% volume CV in the range of 25-55% HCT.

Product storage and handling:

Storage conditions before use: Temperature, 15°C to 28°C; relative humidity, 10% to 90%. Avoid direct sunlight and moisture. Unopened devices may be stored under the above conditions until the expiry date printed on the device packaging.

Post-sampling transport:

It is the responsibility of each analysing laboratory to validate appropriate shipping conditions for the analytes to be tested in that laboratory. It is recommended to use the designated drying pouch for return shipment post-sampling to ensure optimal drying conditions. Our paper-based shipping method has been designed to reduce waste and minimize the carbon footprint of your workflows. However, sensitive analytes may require the use of aluminium-coated plastic bags and the inclusion of desiccant pouches.


Sample preparation and downstream analysis:

Take general safety precautions when handling blood. Samples should be completely dried before removal from the card. Depending on the environment, drying times may vary: drying times between 8 and 24 hours are common.

Remove the peel tabs on the backside of the card to expose the sample discs. Remove the disc using a flat-tip tweezer or a disc removal instrument. When removing the discs manually with tweezers, ensure that most of the paper material is removed: only a thin layer of paper residue may remain on the card.

After collecting the sample discs, proceed with elution. The choice of elution protocol and buffer depends on the target analyte(s) and the downstream analysis platform. The laboratory performing the analysis is responsible for validating assay performance.

Our team of application experts can guide you in the analytical workflow.

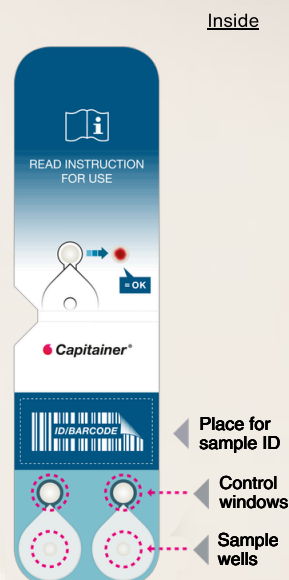

Sampling instructions:

Find Capitainer®B10 Vanadate sampling instructions on our Instructions page [here](#), or scan the QR code in the bottom-left corner.



Capitainer AB, Solna Torg 19, SE-171 45 Solna
info@capitainer.com
www.capitainer.com
COPYRIGHT 2020 © CAPITAINER AB, ALL RIGHTS RESERVED
2026-09-08



www.capitainer.com

Capitainer®