


START  
HERE

# Capitainer® B<sup>50</sup>

Instructions For Use  
Blood Collection Kit

## INCLUDED IN THE SAMPLING KIT.

For single use only.

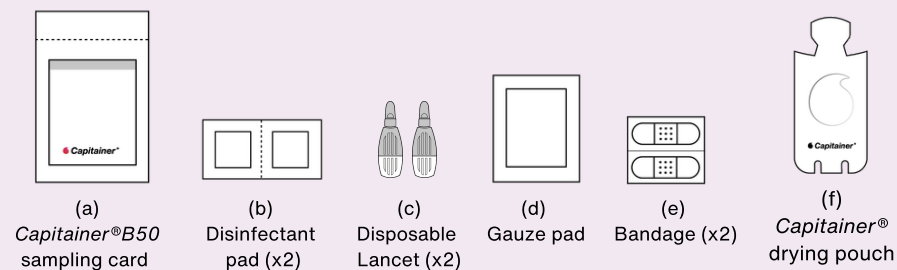

## You will also need:

- A sample ID from your kit provider
- Instructions for how to send your sample to the laboratory



Sample ID

## 01 PREPARE

## 02 COLLECT

## 03 PACK & RETURN

# 01

## PREPARE


1



### Read first

- Read all instructions before you start.
- Scan the QR code to watch the video.


2

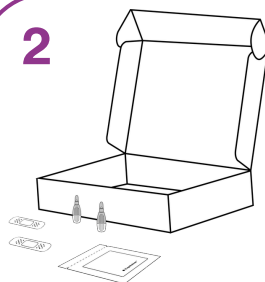

### Prepare material

- Stand or sit at a table.
- Place all items on the table.
- Open the sampling card (a), gauze pad (d), and bandage packages (e).

3

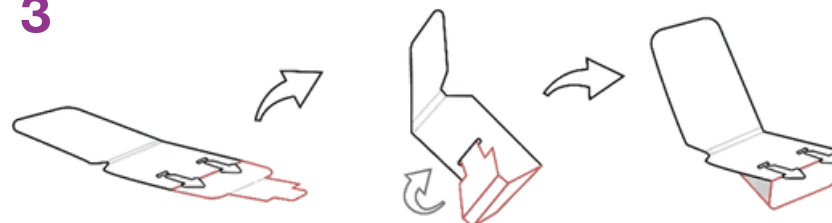

### Prepare card

1

- Open the lid of the sampling card (a). Do not touch the channels.

2

- Fold back the tab to make a stand.

3

- Place the card on the table with the channels facing up.

### Warnings and cautions


- Use the card within 30 minutes of opening.
- Do not bend or twist the card.
- Leave the card on the table while sampling; Do not hold it.

# 02

## COLLECT


02 COLLECT

03 PACK & RETURN

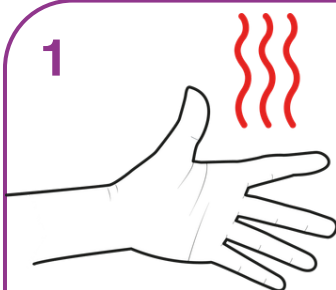


### 1 Warm your hands

- Warm your hand for 2–3 minutes.
- Wash your hands with soap and **warm** water.
- Dry your hands with a clean towel.

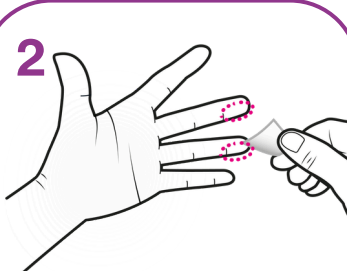


### 2 Prepare finger

- Choose a middle or ring finger.
- Clean the finger with disinfectant pad (b).
- Let the finger dry for 15 seconds.

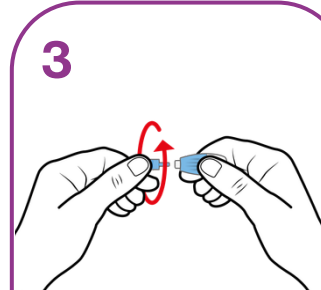

### 3 Prepare lancet

- Twist off the lancet cap (c) (Do Not Pull).

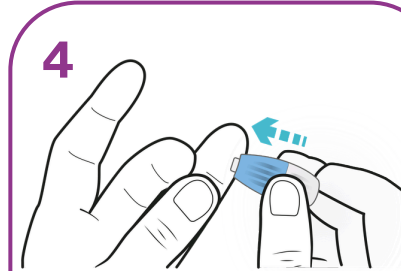


### 4 Prick the finger

- Press your thumb against the fingertip.
- Place the lancet (c) on the side of the fingertip.
- Press firmly to activate the lancet.
- One finger prick should fill both channels.




### 5 Fill one of the channels

- Release thumb pressure and wipe away the first blood drop with the gauze pad (d).
- Let a new **large** blood drop form.
- Fill the channel with blood.
- Make sure the blood goes **into the channel**, not on top.
- Fill to the stop line with as many blood drops as needed.

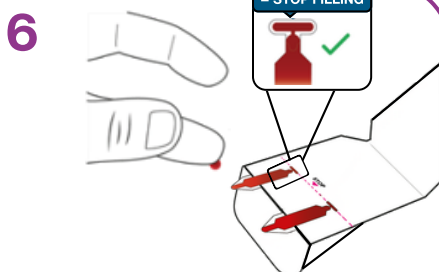

### 6 Fill the other channel

- Fill the other channel to the stop line with the same finger if it is still bleeding. **If not, do a second fingerstick on another finger starting with Step 1.**
- Once both channels have been filled stop the bleeding with the gauze pad (d).
- Apply bandage (e).

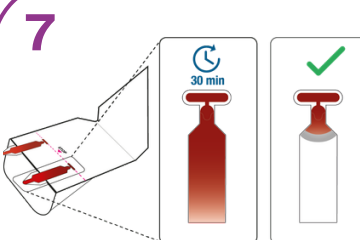


### 7 Check the control window

**Note:** The channels should empty within a few minutes, as shown above. If not, wait 30 minutes before the next step.

# 03

## PACK & RETURN

Throw away all used materials.

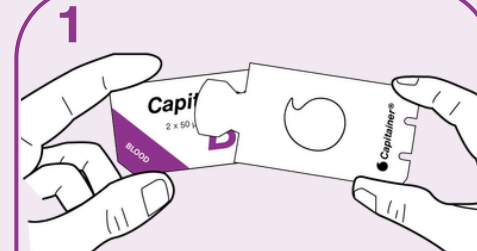

### 1 Finish

- Label the sampling card (a).
- Fold down the stand, remove the sticky peel from the tab and close the protective tab over the channels.
- Place the closed card in the drying pouch (f) and close the flap.
- Return the sample as instructed by your test provider.