


START HERE

Capitainer® DIP70

Instructions For Use
Urine Collection Kit

INCLUDED IN THE SAMPLING KIT.

For single use only.

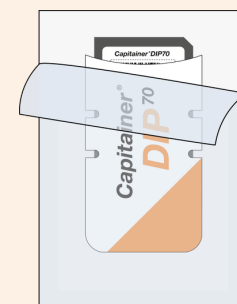


◀ (a) Capitainer®DIP70
sampling card

◀ (b) Capitainer®DIP70
Drying pouch



(c)
Sample
container

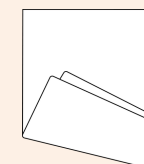

(d)
Absorbent
paper

You will also need:

- A sample ID from your kit provider
- Instructions for how to send your sample to the laboratory



Sample ID

01 PREPARE

02 COLLECT

03 PACK & RETURN

01

PREPARE

1



Read first

- Read all instructions before you start.
- Scan the QR code to watch the video.

2



Wash hands

- Wash your hands with soap and warm water. Dry with a clean towel.
- Choose a place where you can sit or stand.
- Place all the materials on a clean, flat surface.


3



Open packages

- Open the sampling card package.
- Take out the sampling card (a) and the drying pouch (b).

4

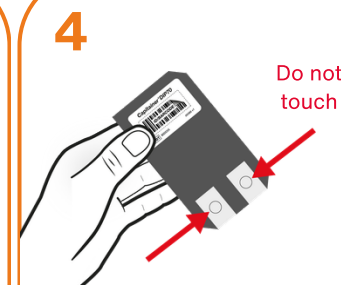

Prepare card

- Place the sampling card (a) on a clean, flat surface without touching the sampling area.

5



Label card

- Label the sampling card (a) with the sample ID.

Warnings and cautions

- Use the card within 30 minutes after opening the package.
- Do not touch the sampling area on the card.

02

COLLECT



Warnings and cautions


- Collect enough urine in the sample container to fully cover the sample discs before dipping.
- Do not go past the cavities when sampling.

02 COLLECT



03 PACK & RETURN


1

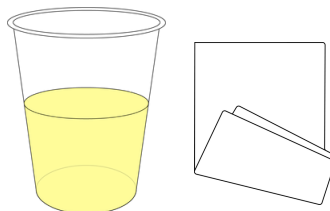

Collect urine

- Have the sample container (c) and the absorbent paper (d) ready.
- Collect urine in the sample container.


2

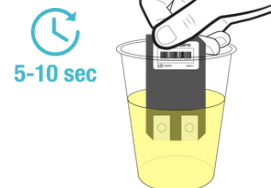

Dip card

- Hold the sampling card (a) upright.
- Dip the sampling discs into the urine for 5–10 seconds.

3

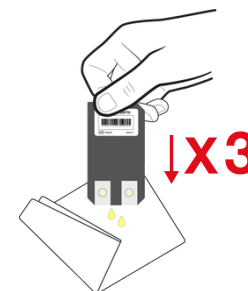

Tap card

- Remove the sampling card (a) from the urine.
- Tap the dipped end firmly on the absorbent paper (d) **three times** to remove extra liquid.


03 PACK & RETURN

03

PACK & RETURN

Throw away all used materials.

1

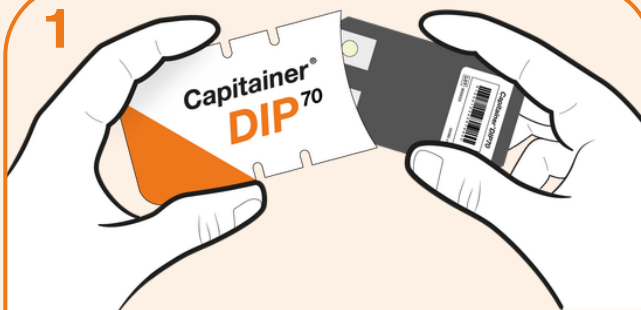

Finish

- Put the card into the drying pouch (b).
- Close the flap.
- Return the sample as instructed by your test provider.