

START
HERE!

SEP¹⁰

For Professional Use Only

Instructions for use.

Sampling kit 2 x 10µl.



INCLUDED IN THE SAMPLING KIT.

Single use products.



(a)
Capitainer® SEP
10 sampling card



(b)
Disinfectant
pad (x2)



(c)
Disposable
lancet (x2)



(d)
Non-adhesive
pad



(e)
Plaster (x2)



(f)
Capitainer®
drying pouch

Additional items you will need:

- A 1-2 cm high box to use as a stand when performing the sampling
- Sample ID

01

PREPARE

Warm hands are critical
for successful sampling.

You can for example:


- Take a shower
- Exercise
- Hold a cup with hot drink for 1-2 mins
- Soak your hand in warm water

1




- Read entire instruction before you start
- Scan QR code or enter <https://capitainer.com/sampling-video-sep10/> to access the video instructions

2

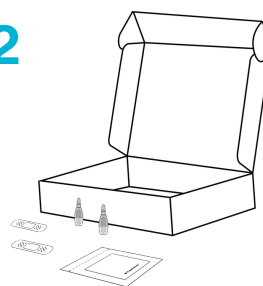


- Choose a place where you can stand comfortably – preferably a dining table. If you prefer to sit, a coffee table works well
- Place all material on table
- Open all packages

3

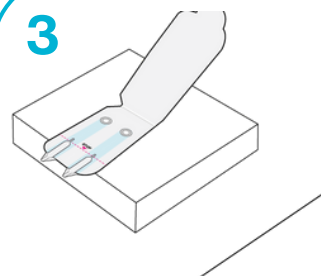

- Open lid of sampling card (a)
- Place card on top of box with sample inlet tips protruding over edge as shown
- Do not touch sample inlets on card

4

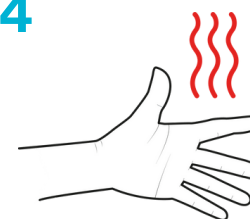

- Choose a hand to sample from and warm for 2-3 mins using a method of choice
- Wash hands with soap and warm water. Dry fully with a clean towel



**Cold hands cause
unsuccessful sampling!**

5



- Swing your arms to increase blood circulation
- During sampling, keep hand positioned below your elbow




Continue on
page two

02


SAMPLE

Warnings and cautions


- Do not touch the sample inlets before sampling
- Use the card within 30 mins of opening the package
- Do not bend or twist the card

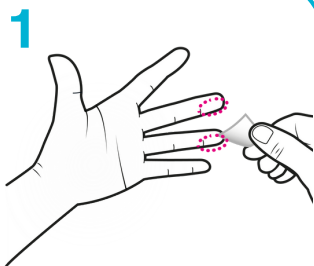


- Choose a finger to prick: middle or ring finger
- Disinfect with disinfectant pad (b) and allow to dry for 15 seconds

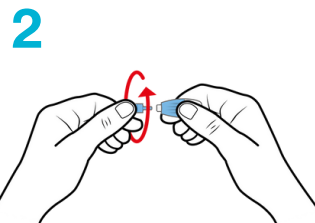


- Twist off and discard lancet (c) cap. DO NOT PULL
- Gently press your thumb against chosen fingertip to fill it with blood

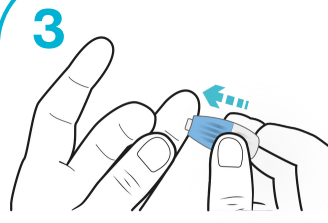

- Place lancet (c) against side of fingertip and press firmly
- Lancet (c) activates automatically and retracts into device
- Release pressure on fingertip and wait for blood to appear

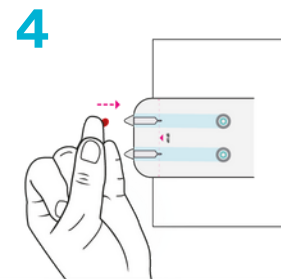


- Wipe off first drop with non-adhesive pad (d)
- Let a **new large drop** of blood form
- Touch the inlet's side with the blood drop, not your finger
- Make sure blood enters into channel and not on top of it

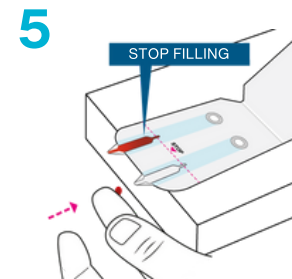


- **Add blood until** it reaches stop line
- Continue to **add** blood into the second channel **until** it reaches stop line
- Do **not** reapply more blood. It may result in an incorrect sample.

6 Important step, read carefully!

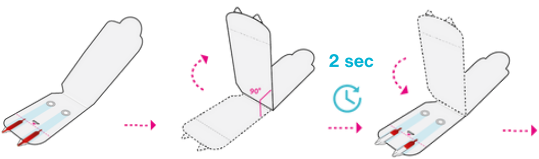

- When sampling is complete: tilt card at 90° angle with inlets facing up for two seconds and put it back flat again.
- Never touch, press or write on the card during or after use; this can ruin your sample

7

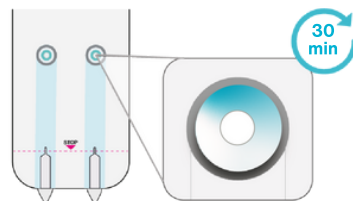


- Leave the card flat, back on the table for 30 minutes.
- Stop bleeding with non-adhesive pad (d)
- Apply plaster (e)
- After 30 minutes, observe that both channels are emptied and that the control windows shows a blue colour.

03

PACK & RETURN

All used materials such as cards, lancets, gauzes, swabs with blood residues must be handled and disposed of safely and in accordance with local regulations.

1

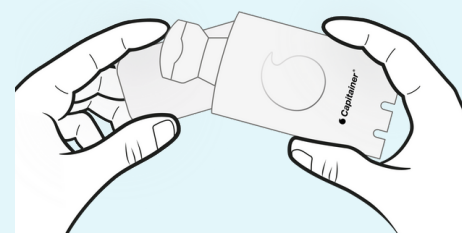

- Label your sampling card
- Place sampling card with lid closed in drying pouch (f)
- Return sample as directed by your test provider.